

lambda250K



Large Area **Medipix3** Based Detector Array

a DESY spinoff company



Designed for high-end X-ray imaging, the LAMBDA 250K provides fast, sensitive detection in a compact package using the Medipix3 chip.

Pioneered by one of the leading research institutes of the world LAMBDA cameras provide the speed and resolution for even the most demanding tasks.




LAMBDA 250K specifications

Sensor	Silicon	GaAs	CdTe
Sensor thickness	300 μm	500 μm	1000 μm
Pixel size	55 μm x 55 μm		
Detector layout (1 module)	1 sensor		
No of Pixels	512 x 512		
Dynamic range	24 bits maximum (dependent on readout mode)		
Energy range	6 – 25 keV	8 – 75 keV	8 – 150 keV
Adjustable threshold range	4 – 40 keV	5 – 50 keV	5 – 75 keV
Frame rate	Up to 2000 Hz (12-bit mode)		
Readout time	No readout time in 12-bit mode, 1 ms in 24-bit mode		
Point spread function	1 pixel FWHM		
External trigger / gate	3.3V TTL		
Software interface	Tango server or open-source hardware library available		
Cooling	Air cooling		
Dimensions (L*H*W)	TBA		
Weight	TBA		
max countrate with correction (10% dev)	2.5×10^8 counts/mm ² /s		

Specifications are subject to change without notice

LAMBDA sensors: Photoelectric absorption

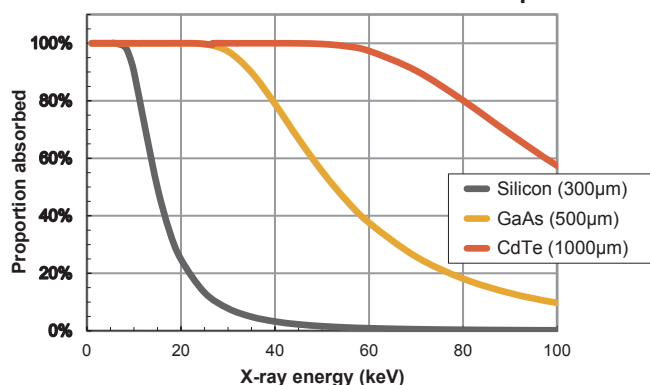

The LAMBDA pixel detector is available with different sensor layers for different X-ray energy ranges. For hard X-ray detection, the GaAs and CdTe LAMBDA systems replace the standard silicon sensor layer in LAMBDA. This provides high quantum efficiency at high X-ray energies (75% at 40keV for GaAs, and 75% at 80keV for CdTe), while retaining all other specifications.

Contact us any time at info@x-spectrum.de or visit us at www.x-spectrum.de

X-Spectrum
enabling research